

 AXIS-SHIELD <small>Innovation for Life</small>	Axis-Shield Active-B12 (Holotranscobalamin) EIA			en
	REF FMABT100		 IVD	

FOR PROFESSIONAL USE ONLY

Read highlighted changes, revised July 2016

INTENDED USE

The Axis-Shield Active-B12 (Holotranscobalamin) assay is an enzyme-immunoassay (EIA) for the quantitative determination of holotranscobalamin (HoloTC) in human serum. HoloTC (vitamin B₁₂ bound to transcobalamin) is used as an aid in the diagnosis and treatment of vitamin B₁₂ deficiency.

INTRODUCTION

Three binding proteins are involved in the transport of vitamin B₁₂ around the body - Intrinsic Factor (IF), transcobalamin (TC) and haptocorrin (HC). These binding proteins ensure the efficient uptake of the very small amounts of vitamin B₁₂ available from the diet. When TC and HC bind vitamin B₁₂ the resulting complexes are known as holotranscobalamin (HoloTC) and holohaptocorrin (HoloHC) to distinguish them from the proteins carrying no vitamin.

The major fraction in the circulation, HoloHC, represents 70-90% of vitamin B₁₂ in the blood but is biologically inert. HoloTC represents only 10-30% of vitamin B₁₂ circulating in the blood but is the only form of vitamin B₁₂ that can be taken up by cells in the body. The TC protein alone transports vitamin B₁₂ from its site of absorption in the ileum to tissues and cells. The vitamin is then internalised as the HoloTC (vitamin B₁₂ bound to transcobalamin) complex via a specific receptor-mediated uptake. This process delivers vitamin B₁₂ into the cells of the body and provides the vitamin as a co-enzyme for essential cellular functions such as DNA synthesis.








As HoloTC has a shorter circulating half-life compared to HoloHC the earliest change that occurs on entering negative vitamin B₁₂ balance is very likely to be a decrease in serum HoloTC concentration¹.

The measurement of Total Serum B₁₂ suffers from some limitations; in particular, most of the measured cobalamin is bound to biologically inert HC. Several studies have been published which conclude that HoloTC would be a better indicator of vitamin B₁₂ status than Total Serum B₁₂^{2,3}. As expected, HoloTC levels are low in patients with biochemical signs of vitamin B₁₂ deficiency⁴. Low values have been reported in vegetarians⁵, vegans⁶, and in populations with a low intake of vitamin B₁₂⁷. Notably, low levels of HoloTC but not Total B₁₂ in serum were reported in patients with Alzheimer's disease compared to levels in a healthy control group⁸. HoloTC levels reflect vitamin B₁₂ status, independent of recent absorption of the vitamin⁹.

PRINCIPLE OF THE ASSAY

The microtitre wells are coated with a highly specific monoclonal antibody for Active-B12 (Holotranscobalamin). During the first incubation holotranscobalamin in serum specifically binds to the antibody-coated surface. In the second incubation the Conjugate binds to any captured holotranscobalamin. The wells are then washed to remove unbound components. Bound holotranscobalamin is detected by incubation with the Substrate. Addition of Stop Solution terminates the reaction, resulting in a coloured end-product. The concentration of holotranscobalamin in pmol/L is directly related to the colour generated and can be estimated by interpolation from a dose-response curve based on Calibrators.

KIT COMPONENTS

CONJ	1 × 15 mL	Alkaline phosphatase-labelled murine monoclonal antibody to human transcobalamin in Tris buffer with protein stabiliser. Preservative: < 0.1% (w/v) sodium azide. Ready-to-use.	
SUBS	1 × 15 mL	para-NitroPhenyl Phosphate (pNPP), buffer solution. Ready-to-use. Do not expose to light during storage. N.B. WARNING.	 
SOLN STOP	1 × 15 mL	1M Sodium hydroxide, (pH >10). Ready-to-use. N.B. DANGER.	 
BUF WASH 8X	2 × 25 mL	Phosphate buffer. Preservative: 0.72% (w/v) sodium azide. Dilute before use. N.B. WARNING.	
MTP 8X12	12 × 8 well microtitre (breakapart) strips	Coated with anti-holotranscobalamin murine monoclonal antibody, in a resealable foil pack with desiccant.	
CAL A - CAL F	6 × 1.0 mL	Cal A is a phosphate buffer with protein (bovine) stabiliser. Cals B-F are phosphate buffer with protein (bovine) stabiliser containing HoloTC. Preservative: < 0.1% (w/v) sodium azide. Ready-to-use. Do not expose to light during storage. SEE VIAL LABELS FOR CONCENTRATIONS.	
CONTROL L	1 × 1.0 mL	Phosphate buffer with protein (bovine) stabiliser containing HoloTC. Preservative: < 0.1% (w/v) sodium azide.	
CONTROL H	1 × 1.0 mL	Ready-to-use. Do not expose to light during storage.	
PRE-TREATMENT	1 × 25 mL	Citrate buffer. Preservative: < 0.1% (w/v) sodium azide. Ready-to-use.	

STANDARDISATION

There is currently no internationally recognized reference method or reference material for standardization. The Axis-Shield Active-B12 (Holotranscobalamin) Calibrators are traceable to internal reference standards which underwent a one-time value assignment.

STORAGE OF REAGENTS

Opened (In-Use) Kit Stability

A kit was opened, and reused on three occasions over a three month period with no adverse effect on kit performance. Following use, components must be returned to storage at 2-8°C.

Unopened kit stability

At 2-8°C unopened components are stable until as directed on the labels.

Handling and Procedural Notes

1. Store kit components at 2-8°C and use until the expiry date on the labels. Do not use expired reagents.
2. Each lot of reagents and calibrators has been standardised to produce the correct reaction. Do not interchange the reagents or calibrators between lots.
3. Calibrator concentrations are displayed on vial labels and may vary between lots.
4. Do not freeze kits.
5. Wash Buffer Concentrate must be diluted before use. All other reagents are ready-to-use.
6. Diluted Wash Buffer is stable for at least 3 months if microbial contamination is avoided. Return to 2-8°C storage after each use.
7. Replace surplus (unused) microtitre strips in the foil pack with the desiccant. Ensure seal is integral and return to 2-8°C, until required.
8. Do not expose the Calibrators, Controls or Substrate to light during storage.
9. Avoid contamination of reagents. Use a new disposable pipette tip for each reagent or sample manipulation.

Indications of Deterioration

The Substrate should be colourless to pale yellow in colour. Darker yellow colouring indicates contamination and the reagent must be discarded. Turbidity or precipitation in any component indicates deterioration and the component should be discarded.

Sample Collection and Storage



1. The assay is recommended for human serum (including serum separator tubes).
2. Do not use grossly haemolysed or turbid samples.
3. Thoroughly mix thawed samples before assay and avoid repeated freeze/thawing.
4. Samples may undergo 3 freeze thaw cycles. Thawed samples should be centrifuged at $\geq 10,000g$ for 5 minutes before assaying.
5. Do not subject on-the clot or off-the clot samples to temperature above room temperature for longer than overnight (≤ 16 hours).
6. Samples may be stored at 2-8°C on-the clot for up to 3 days or off-the-clot for four weeks; for longer storage samples must be stored off-the-clot at -20°C for up to 6 months.
7. Prepare each sample prior to assay by adding an equal volume of Pre-Treatment to Sample e.g. 150 μ L sample plus 150 μ L Pre-Treatment. Pre-treated samples may be stored capped for up to 24 hours at 2-8°C prior to assay.

WARNINGS AND PRECAUTIONS

For *in vitro* diagnostic use only. Safety Precautions

1. Adhere strictly to the instructions in this booklet, particularly for handling and storage conditions.
2. Do not pipette by mouth.
3. Do not smoke, eat, drink or apply cosmetics in areas where kits and samples are handled.
4. Any skin complaints, cuts, abrasions and other skin lesions should be suitably protected.
5. The Calibrators, Controls, Conjugate, Pre-Treatment and Wash Buffer Concentrate contain sodium azide which can react with lead and copper plumbing to form highly explosive metal azides. On disposal, drain with large quantities of water to prevent azide build-up.
6. The Stop Solution contains sodium hydroxide. Avoid contact with skin, eyes and mucous membranes. Spillage should be mopped up with copious amounts of water. If contact with skin or eyes occurs, irrigate with water and seek medical attention immediately.
7. Material safety data sheets for all components contained in this kit are available on request from Axis-Shield Diagnostics.
8. This product requires the handling of specimens and materials of human and animal origin. It is recommended that all human and animal sourced materials be considered potentially infectious and be handled in accordance with the OSHA Standard on Bloodborne Pathogens Biosafety Level 2 or other appropriate biosafety practices should be used for materials that contain or are suspected of containing infectious agents.

Caution: Federal law restricts this device to sale by or on the order of a physician.

 Danger	Stop Solution	<p>WARNING H314 – Causes severe skin burns and eye damage.</p> <p>PREVENTION P264 – Wash hands thoroughly after handling. P280 – Wear protective gloves/protective clothing/eye protection/face protection.</p> <p>RESPONSE P305+351+338 – IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.</p>	
 Warning	Substrate	<p>WARNING H302 – Harmful if swallowed. H312 – Harmful in contact with skin. H332 – Harmful if inhaled.</p> <p>PREVENTION P260 – Do not breathe dust/fume/gas/mist/vapours/spray. P271 – Use only outdoors or in a well-ventilated area. P280 – Wear protective gloves/protective clothing/eye protection/face protection.</p> <p>RESPONSE P301+310 – IF SWALLOWED : Immediately call a POISON CENTRE or doctor/physician. P304+340 – IF INHALED : Remove victim to fresh air and keep at rest in a position comfortable for breathing. P312 – Call a POISON CENTRE or doctor/physician if you feel unwell.</p>	
	Wash Buffer Concentrate (8X)	<p>WARNING H302 – Harmful if swallowed. H412 – Harmful to aquatic life with long lasting effects. EUH032 – Contact with acids liberates very toxic gases.</p> <p>PREVENTION P264 – Wash hands thoroughly before handling. P270 – Do not eat, drink or smoke when using this product. P273 – Avoid release to the environment.</p> <p>RESPONSE P301+310 – IF SWALLOWED : Immediately call a POISON CENTRE or doctor/physician. P330 – Rinse mouth.</p>	

PREPARATION

Materials/Equipment Required but not Provided

1. 96 well plate/strip reader with 405nm filter.
2. Precision pipette(s) to dispense 100 µL. An 8-channel dispenser, or similar, to dispense approximately 250-300 µL for manual washing (for example the StatMatic 8-channel dispenser from Tricontinent).
3. Glass/plastic measuring cylinder 1x200 mL.
4. Distilled/deionised water.
5. Paper towels.
6. Timer for 30, 35 and 60 minute intervals.

Preparation for the Assay

Allow all kit components, including the microtitre strips, to warm up to 18-25°C for 30-60 minutes before use. Mix reagents by gentle inversion.

When stored at 2-8°C the wash buffer will precipitate (crystals may be visible). Before diluting in water, allow the wash buffer to warm up (can be placed in an incubator at 37°C if required to speed the process) until **NO** precipitation is evident to the naked eye.

Dilute the following reagent and mix thoroughly:

Reagent	Volume	Add
Wash Buffer Concentrate x 8	1 vial	175 mL distilled/deionised water

Calculate the number of microtitre strips required for the current assay and retain these in the microtitre strip holder. Return surplus strips to the resealable foil pack with the desiccant and store at 2-8°C until required. Ensure that all strips are securely held within the microtitre strip holder. Users may wish to number each strip along the top edge to aid identification. Retain the microtitre strip holder for future use.

Prepare each sample prior to assay by adding an equal volume of Pre-Treatment to sample e.g. 150 µL sample plus 150µL Pre-Treatment. Pre-treated samples may be stored capped for up to 24 hours at 2-8°C prior to assay.

ASSAY PROTOCOL


1. Reference wells for identification.
2. Pipette 100 μ L Calibrators in duplicate, Kit Controls in duplicate and pre-treated (50:50) patient samples in duplicate, into appropriate wells. Remember to change pipette tips between additions. This step should not exceed 15 minutes.
3. Incubate 60 ± 10 minutes at $18-25^{\circ}\text{C}$.
4. Decant strip contents by quick inversion over a sink suitable for the disposal of biological materials, bearing in mind the potential infective hazard of the samples. Blot inverted strips well with paper towels. **Do not wash.**
5. Add 100 μ L Conjugate to each well.
6. Incubate 35 ± 5 minutes at $18-25^{\circ}\text{C}$.
7. Decant strip contents by quick inversion over a sink suitable for the disposal of biological materials. Blot inverted strips well with paper towels.
8. Wash wells **five times** with a minimum of 250 μ L diluted Wash Buffer. **Decant and blot after each wash addition.**
9. Add 100 μ L Substrate to each well.
10. Incubate 30 ± 5 minutes at $18-25^{\circ}\text{C}$. **Do not decant.**
11. Add 100 μ L Stop Solution to each well, in the same order and rate as the Substrate. Tap wells gently to mix.
12. Read strips at 405nm. Read within 120 minutes of addition of Stop Solution.

CALCULATION and interpretation OF RESULTS

Plot the mean absorbance value of each Calibrator on the y-axis against the corresponding concentration in pmol/L on the x-axis. Calibrator concentrations are displayed on vial labels. CONCENTRATION VALUES ARE ASSIGNED TO EACH LOT OF CALIBRATORS and may vary between lots.

The concentration (pmol/L) of each sample can be calculated by locating the point on the curve corresponding to the mean sample absorbance value and reading the corresponding concentration in pmol/L from the x-axis. This procedure can be performed manually using graph paper or using a plate reader with software incorporating curve fitting procedures. If using a plate reader with internal software, a linear regression curve-fit algorithm should be used.

A typical calibration plot is shown below for reference purposes, it must not be used for interpreting results. Samples with concentrations above 128 pmol/L are outside the range of the assay, and should be recorded as $>128\text{pmol/L}$ and results must not be extrapolated. Individual sample replicates deviating less than 20% can be taken to indicate acceptability of the assay.

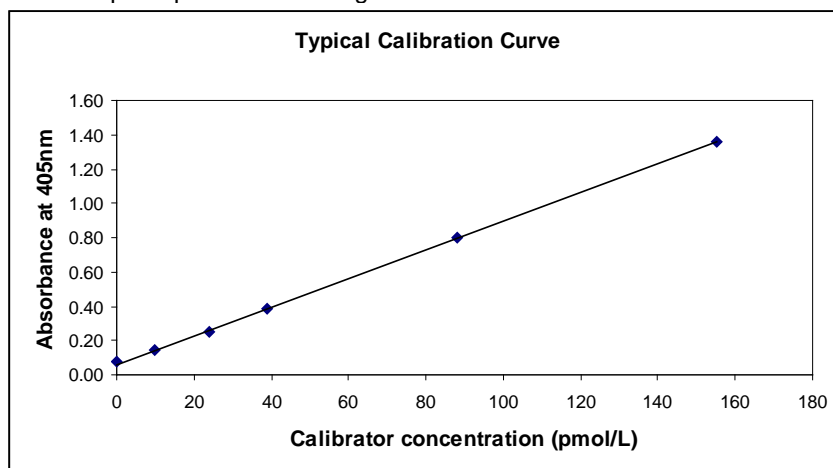

RESULTS

Unit of Measure

The unit of measurement for the Axis-Shield Active-B12 (Holotranscobalamin) assay is pmol/L.

Measuring Interval (Reportable Range)

The measurable range of the assay is 10 pmol/L to 128 pmol/L.

QUALITY CONTROL

Ensure that adequate maintenance and calibration of the plate-reader is performed according to the manufacturer's instructions, and that the correct wavelength (405nm) and curve-fit algorithm (linear regression) are employed.

Users should ensure they are fully acquainted with the instructions for the assay, particularly the Warnings and Precautions section, and the Handling and Procedural Notes. Users should demonstrate that they can obtain performance specifications for precision and reportable range of test results comparable to those established by the manufacturer before reporting patient test results. Ready-to-use Low and High Kit Controls must be run in duplicate in all assays to monitor the quality of the test procedure.

Assuming the precision specifications described by the manufacturer are met, failure of any Control to meet the Control specifications below renders the assay invalid and patient results should not be reported. The operator may repeat the assay, having reviewed their procedure, or contact the manufacturer. If repeating the assay, prepare a fresh dilution of each sample. Laboratories may wish to include in-house controls in each assay run. Store such control material at or below -20°C and avoid repeat freeze/thaw cycles. Preservatives such as sodium azide at <0.1% (w/v) will not affect sample results. Reference ranges and appropriate cut-off points should be calculated for the specific populations served by users.

Control Specifications

	Low	High
Specification (mean of duplicate)	15 to 35 pmol/L	36 to 84 pmol/L

EXPECTED VALUES

135 serum samples from asymptomatic apparently healthy donors with an age range of 18-75 years, comprising approximately equal numbers of males [n = 65] and females [n = 70], were tested in the Axis-Shield Active-B12 (Holotranscobalamin) EIA. The overall mean Active-B12 (Holotranscobalamin) concentration for this population was 72 pmol/L (range 15 to 147 pmol/L). On the basis of this reference population data, the reference range (central 95% of the results) is:

Reference Range

21 – 123 pmol/L

This reference range is suggested as a guideline only and each laboratory should establish its own reference range which may be unique to the population it serves depending upon geographical, patient, dietary or environmental factors or clinical practice.

PERFORMANCE DATA

Representative data; results in individual laboratories may vary.

Dilution Linearity

Based guidance from CLSI document EP6-A10, the Axis-Shield Active-B12 (Holotranscobalamin) EIA demonstrated linearity across the measuring range of the assay as demonstrated in a study from 5.3 to 156.0 pmol/L (rounded to 1 decimal place).

Accuracy

A correlation study was performed with serum specimens from apparently healthy adults. All specimens were analysed using the Axis-Shield Active-B12 (Holotranscobalamin) EIA and another commercially available Holotranscobalamin assay according to the CLSI document EP9-A211. Specimen concentrations ranged from 13.8 to 112.8 pmol/L in the assay. The data obtained gave the following statistical values:

Axis-Shield Active-B12 (Holotranscobalamin) EIA versus a commercially available assay	
Number of specimens	111
Slope of regression line (Passing-Bablok regression) (95% CI)	0.95 (0.89 to 1.01)
Y-intercept (Passing-Bablok regression) (95% CI)	8.39 (5.73 to 11.77)
Correlation coefficient (r) (Pearson) (95% CI)	0.93 (0.90 to 0.95)

Precision

6 human serum samples were assayed at Axis-Shield using 3 lots of reagents. Samples were assayed by 2 operators in replicates of 8, once a day for 5 days (total n=80). Data from this study are summarised in the following table:

Sample	n	Lot	Operator	Mean (pmol/L)	Intra-assay %CV	Total %CV
1 A	80	1	1	17.8	7.5%	8.2%
			2	17.5	3.1%	9.3%
		2	1	20.1	6.0%	6.6%
			2	20.3	6.9%	9.2%
		3	1	19.1	5.5%	8.0%
			2	18.9	8.5%	11.0%
2 A	80	1	1	21.8	5.5%	9.9%
			2	21.8	3.9%	7.5%
		2	1	22.6	5.6%	8.7%
			2	23.5	9.0%	10.3%
		3	1	23.9	7.0%	10.2%
			2	23.2	5.8%	8.9%
3 A	80	1	1	28.8	3.8%	7.8%
			2	30.7	4.3%	9.6%
		2	1	31.0	6.8%	8.0%
			2	31.4	4.3%	6.1%
		3	1	31.5	4.5%	6.4%
			2	32.2	4.0%	9.2%
4 A	80	1	1	49.3	3.9%	7.4%
			2	52.6	4.1%	6.7%
		2	1	50.8	5.6%	10.0%
			2	51.7	4.7%	5.9%
		3	1	52.6	4.6%	4.8%
			2	55.0	5.5%	6.1%

5 A	80	1	1	68.4	4.0%	7.6%
			2	73.2	3.7%	7.5%
		2	1	74.8	4.3%	8.2%
			2	75.9	4.6%	6.4%
		3	1	75.1	4.4%	7.9%
			2	76.3	4.9%	6.2%
7 A	80	1	1	115.9	4.2%	5.9%
			2	121.1	3.6%	7.0%
		2	1	123.2	4.3%	10.2%
			2	124.0	4.2%	6.4%
		3	1	127.0	4.8%	10.1%
			2	129.5	3.2%	5.6%
Low Control	80	1	1	23.7	9.4%	10.9%
			2	23.8	5.1%	11.5%
		2	1	20.0	6.0%	7.5%
			2	18.6	5.8%	8.5%
		3	1	20.3	8.3%	9.7%
			2	20.1	8.3%	10.0%
High Control	80	1	1	61.2	6.3%	6.4%
			2	58.8	4.5%	8.9%
		2	1	50.3	6.3%	8.1%
			2	50.2	5.9%	8.4%
		3	1	52.2	7.7%	9.2%
			2	50.8	5.8%	8.5%

Limit of Blank

In a representative study, Limit of Blank determinations were performed using two low-level holotranscobalamin samples and two reagent lots (120 replicates per reagent lot). The limit of blank of the Axis-Shield Active-B12 (Holotranscobalamin) EIA was found to be 4.9 pmol/L (rounded to 1 decimal place).

Limit of Detection

In a representative study, Limit of Detection determinations were performed using five low-level holotranscobalamin samples and two reagent lots (120 replicates per reagent lot). The limit of blank of the Axis-Shield Active-B12 (Holotranscobalamin) EIA was found to be 8.1 pmol/L (rounded up to 1 decimal place).

Limit of Quantitation

Limit of quantitation determinations were performed using five low-level holotranscobalamin samples and two reagent lots (120 replicates per reagent lot). The limit of quantification of the Axis-Shield Active-B12 (Holotranscobalamin) EIA was found to be 8.3 pmol/L (rounded to 1 decimal place).

High Dose Hook

High dose hook is a phenomenon whereby very high level specimens may read within the dynamic range of the assay. For the Axis-Shield Active-B12 (Holotranscobalamin) EIA. No high dose hook effect was detected with two samples with a concentration of approximately 419 and 2236 pmol/L.

Cross-reactivity

The Axis-Shield Active-B12 (Holotranscobalamin) EIA is designed to have a maximum deviation in holotranscobalamin concentration of $\leq 10\%$ in the presence of apotranscobalamin or haptocorrin.

A study was performed based on guidance from the Clinical and Laboratory Standards Institute (CLSI) document EP7-A2¹². Three samples with holotranscobalamin levels across the assay range were supplemented with 500 pmol/L apotranscobalamin or 5000 pmol/L haptocorrin. The maximum deviation in holotranscobalamin concentration ranged from -5% to 1%.

Interference

The Axis-Shield Active-B12 (Holotranscobalamin) EIA is designed to have a maximum deviation in holotranscobalamin concentration of $\leq 10\%$ in the presence of potentially interfering compounds.










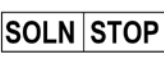





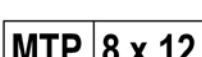

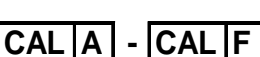


A study was performed based on guidance from the Clinical and Laboratory Standards Institute (CLSI) document EP7-A2¹². Samples with holotranscobalamin levels across the assay range were supplemented with the potentially interfering compounds listed in the table below. The maximum deviation in holotranscobalamin concentration ranged from -10% to 8%.

Potential Interfering Substance	No interference found up to the following concentration
Haemoglobin	500 mg/dL
Bilirubin	30 mg/dL
Triglyceride (Intralipid Solution)	3000 mg/dL
Rheumatoid Factor	7500 IU/dL
Total Protein	9000 mg/dL

REFERENCES

1. Nexø E, Hvas A-M, Bleie Ø *et al.* Holo-transcobalamin is an early marker of changes in cobalamin homeostasis. A randomized placebo-controlled study. *Clin Chem* 2002;48(10):1768-71
2. Valente E, Scott JM, Ueland PM *et al.* Diagnostic accuracy of holotranscobalamin, methylmalonic acid, serum cobalamin, and other indicators of tissue vitamin B12 status in the elderly. *Clin Chem* 2011;57(6):856-863
3. Nexø E, Hoffmann-Lucke E. Holotranscobalamin, a marker of vitamin B12 status: analytical aspects and clinical utility. *Am J Clin Nutr* 2011;94(1):359S-365S
4. Obeid R, Jouma M, Hermann W. Cobalamin status (holotranscobalamin, methylmalonic acid) and folate as determinants of homocysteine concentration. *Clin Chem* 2002;48(11):2064-5
5. Herrmann W, Schorr H, Obeid R *et al.* Vitamin B12 status, particularly holotranscobalamin II and methylmalonic acid concentrations and hyperhomocysteinemia in vegetarians. *Am J Clin Nutr* 2003;78:131-6
6. Lloyd-Wright Z, Hvas A-M, Møller J *et al.* Holotranscobalamin as an indicator of dietary vitamin B12 deficiency. *Clin Chem* 2003;49(12):2076-8.
7. Refsum H, Yajnik CS, Gadkari M *et al.* Hyperhomocysteinemia and elevated methylmalonic acid indicate a high prevalence of cobalamin deficiency in Asian Indians. *Am J Clin Nutr* 2001;74:233-41
8. Refsum H, Smith AD. Low Vitamin B12 status in confirmed Alzheimer's disease as revealed by serum holotranscobalamin. *J Neurol Neurosurg Psychiatry* 2003;74:959-61
9. Chen X, Remacha AF, Sarda MP *et al.* Influence of cobalamin deficiency compared with that of cobalamin absorption on serum holo-transcobalamin II. *Am J Clin Nutr* 2005;81:110-14
10. The National Committee for Clinical Laboratory Standards (NCCLS). *Evaluation of the Linearity of Quantitative Measurement Procedures: A Statistical Approach; Approved Guideline*. NCCLS document EP6-A. Wayne, PA: NCCLS; 2003.
11. The National Committee for Clinical Laboratory Standards (NCCLS). *Method Comparison and Bias Estimation Using Patient Samples; Approved Guideline—Second Edition*. NCCLS Document EP9-A2. Wayne, PA: NCCLS; 2002.
12. Clinical and Laboratory Standards Institute (CLSI). *Interference Testing in Clinical Chemistry; Approved Guideline—Second Edition*. CLSI document EP7-A2. Wayne, PA: CLSI; 2005.

KEY TO SYMBOLS

	<i>In vitro</i> diagnostic medical device		Low Control
	Catalogue number		High Control
	Lot		Conjugate
	96 tests		Substrate
	Caution		Stop Solution
	Consult instructions for use		Pre-Treatment
	Protect from Light		Wash Buffer
	Use by		Microtitre (Breakapart) Strips
	Store at 2-8°C		Calibrators A to F
	For Prescription Use Only		Manufactured by

Axis-Shield Diagnostics Ltd.,
The Technology Park,
Dundee,
DD2 1XA, UK

Tel: +44 (0) 1382 422000
Fax: +44 (0) 1382 422088
e-mail: axd.axis-shield@alere.com
Web: www.axis-shield.com

Ver: 2016/07